


Carnation Ambulatory Monitor (CAM[™])



Actual device size is 7" (long), 1 1/2" (at widest point) 3/8" (deep) and weighs 18 grams.

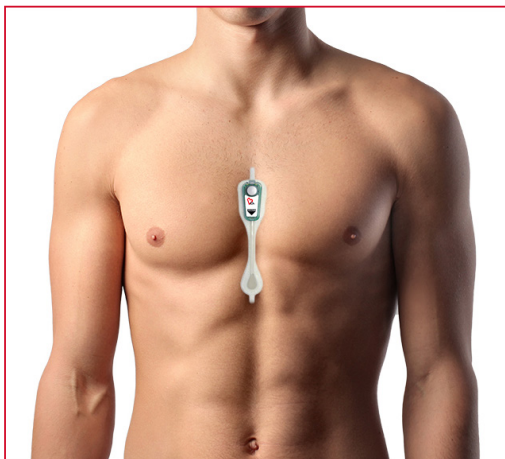

Discreet



Water Resistant



No Wires

The CAM device is bandage size, lightweight, and easy to use. It can be worn during exercise, sleep, or while showering.

CAM is comfortable for extended wear for men and women and has a unique, narrow, flexible, hourglass shape that enables full physical movement.

We're close to your heart[™]





CAM™ Patient Instructions

- Do not shower or engage in any physical activity that will cause sweating for at least 2 hours after CAM is placed on the skin.
- While showering, do not stand with the water stream directly on the CAM.
- DO NOT bathe in a tub, jacuzzi or swimming pool where the CAM is completely under water. Do not use a steam room/sauna or do “hot” yoga.
- Use caution when fastening/unfastening seat belts so as not to rub or dislodge the monitor.
- You are able to exercise with the CAM, however, excessive sweating may cause the CAM to come off early – If so, simply dry the chest area and replace the CAM.
- Please write down BOTH date and time of the event in the diary.
- Press the event button on the CAM if you feel any symptoms.
- Press button once, briefly, each time symptoms occur. Do not fondle or press repeatedly for the same event.

TO REMOVE

- Pull on either tab located at the top or bottom of the CAM – Use adhesive removal, if necessary, to remove extra adhesive – soap and water can be used as a substitute.
- Stick the monitor on the box insert, then place CAM monitor and diary inside the return box.
- Return the monitor as instructed by your physician’s office.

For Any Questions During Your CAM Wear
Please Contact Your Physician

OR

BardyDx Customer Service: **+1-844-422-7393**